


Mercury Medical[®]

Announcing
The New Disposable

air-Q[®] **Blocker**

Masked Laryngeal Airway

*The new design, with integrated channel to access the esophagus, has evolved into a superior airway. Clinicians can now place an NG Tube through the air-Q **Blocker** channel.*

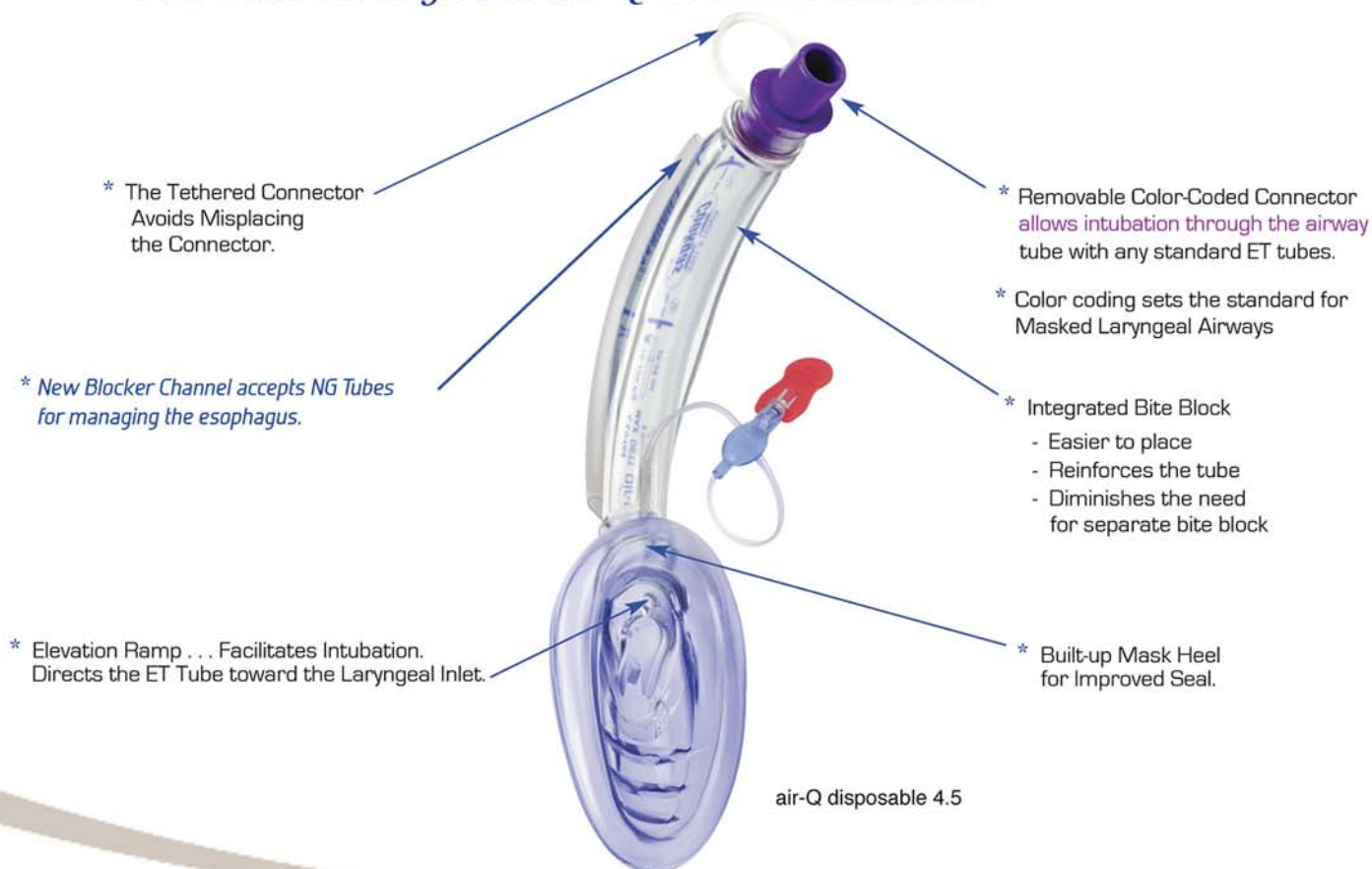

11300 - 49th Street North
Clearwater, Florida USA 33762-1807
www.mercurymed.com

Mercury Medical[®]

New, Disposable air-Q[®] Blocker

Masked Laryngeal Airway

*Developed and Designed
by an Anesthesiologist . . .*

for Use with Confidence!

- *Now with an integrated channel for use with an NG Tube for managing potential gastric fluids.*



The new and improved design has evolved into a superior airway providing:

- *Available in 3 Sizes*
- *Ease of Placement*
- *Better Seal in “Everyday” Use*
- *Improved Intubation Success*

| Product Advantages | | air-Q | Other Airways |
|---|--|-------|---------------|
| <i>New Blocker Channel</i> | <i>Allows for use with an NG tube for managing the esophagus</i> | Yes | No |
| Heel of mask cuff | Built-up for improved seal | Yes | No |
| Elevation Ramp (built into airway outlet.) | Facilitates placement, directs ET tube toward laryngeal inlet. | Yes | No |
| Tethered connector to tube | To avoid misplacing the connector while in use. | Yes | No |
| Built-In Bite Block | -Easier to place - Reinforces the tube - Diminishes the need for separate bite block | Yes | Some |
| Mask tip with recessed front that approximates shape of posterior pharynx. | Easier positioning and improved airway alignment. | Yes | No |
| Specially-designed ridges inside mask cavity. | Provides lateral stability to help prevent tip from bending backward and improves mask seal. | Yes | No |
| Keyhole-shaped airway outlet. | Elevates epiglottis to provide unobstructed access to laryngeal inlet. | Yes | No |
| Removable, tethered connector for direct access to airway tube. | Allows intubation through tube with any standard ET tube. | Yes | No |
| Facilitates intubation, eliminates duplicate inventories and costly special ET tubes. | Cost-effective. | Yes | Some |

Distributed by:

Mercury Medical[®] ANESTHESIA

Clearwater, FL 33762-4807 · 800.237.6418 · www.mercurymed.com

3 sizes available



2.5 disposable
air-Q Blocker



3.5 disposable
air-Q Blocker



4.5 disposable
air-Q Blocker

Ordering Information

| Item # | Size | Description | Quantity | Patient Weight |
|----------|------|---|----------|----------------|
| #10-5025 | 2.5 | air-Q Blocker Disposable Airway only, yellow connector | 10/Box | 30 - 50 kgs |
| #10-5035 | 3.5 | air-Q Blocker Disposable Airway only, red connector | 10/Box | 50 - 70 kgs |
| #10-5045 | 4.5 | air-Q Blocker Disposable Airway only, purple connector | 10/Box | 70 - 100 kgs |
| #10-1005 | #0 | Disposable Removal Stylet for 2.0 & 2.5 airway | 10/Box | |
| #10-1004 | #1 | Disposable Removal Stylet for 3.5 & 4.5 airway | 10/Box | |
| #10-5026 | 2.5 | air-Q Blocker Disposable Airway kit w/lubricant & 10 cc syringe with swivel elbow, yellow connector | 10/Box | 30 - 50 kgs |
| #10-5036 | 3.5 | air-Q Blocker Disposable Airway kit w/lubricant & 10 cc syringe with swivel elbow, red connector | 10/Box | 50 - 70 kgs |
| #10-5046 | 4.5 | air-Q Blocker Disposable Airway kit w/lubricant & 10 cc syringe with swivel elbow, purple connector | 10/Box | 70 - 100 kgs |

Registered trademark, air-Q®, property of Cookgas, LLC

Patents U.S.

Patents Canada

Patents UK

5,937,860
US 6,422,239 B1
US 6,892,731 B2
US 7,331,347 B2
US 6,705,321 B2
US 7,357,845 B2
US 7,780,900 B2

2,231,331
2,371,435

GB2357437B
GB2407293B
GB2405589B
GB2324040B

Other patents pending.



... The everyday airway designed to facilitate intubation.

PRSRT STD
U.S. POSTAGE
PAID
Permit #579
St. Petersburg, FL

MercuryMedical[®]

air-Q[®] *Blocker*

*... The everyday disposable airway
designed to facilitate intubation.*

air-Q **Blocker** Disposable 3.5 Masked Laryngeal Airway
Small Adult with an integrated channel to the
esophagus.

Single-Patient-Use

