




DIMAR

Medical Devices



DimAir Line

**Ambulance
Emergency
Department**

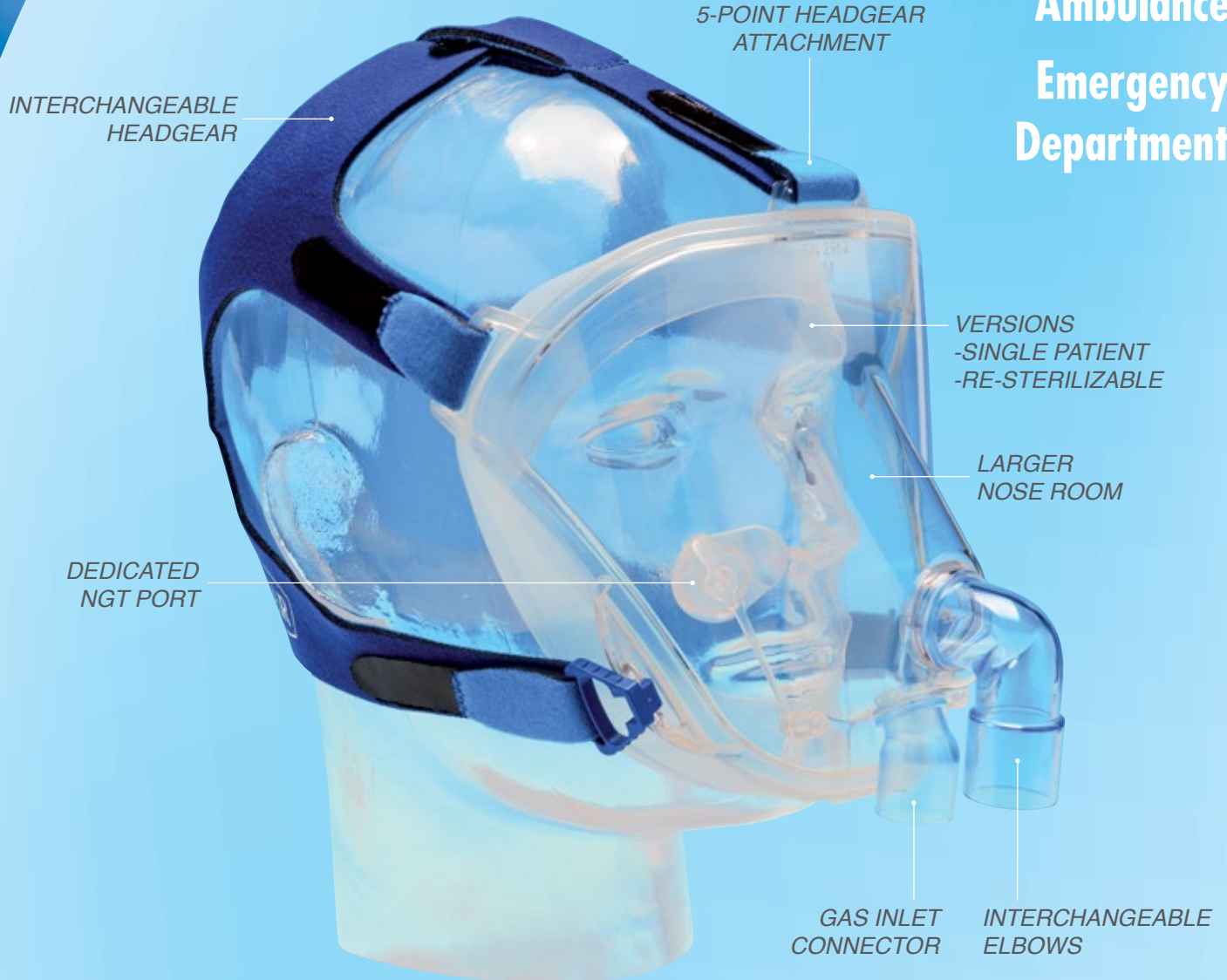

INTERCHANGEABLE HEADGEAR

5-POINT HEADGEAR ATTACHMENT

VERSIONS
-SINGLE PATIENT
-RE-STERILIZABLE

LARGER NOSE ROOM

DEDICATED NGT PORT

GAS INLET CONNECTOR

INTERCHANGEABLE ELBOWS

TOTAL FACE DIMAX ZERO


“DiMax Zero” TOTAL FACE MASK

A new generation total face mask for Non-Invasive Ventilation in PSV, BIPAP and CPAP mode. Its original configuration reduces to zero the CO₂ rebreathing. Available in different sizes to fit any paediatric or adult patient, DiMax can also be used in combination with dual or single limb ventilators as well as CPAP high flow generators by simply changing the elbow. Unique features such as the 5-point mask attachment, NGT airtight port, bronchoscopy airtight port, extra nose room, make DiMax Zero the perfect interface for Non Invasive Ventilation.

8 plus

- zero rebreathing
- securement points
- comfort
- sealing
- practicality
- versatility
- safety
- versions

Fluid Dynamics of CO₂ rebreathing in standard total face mask compared to DiMax Zero

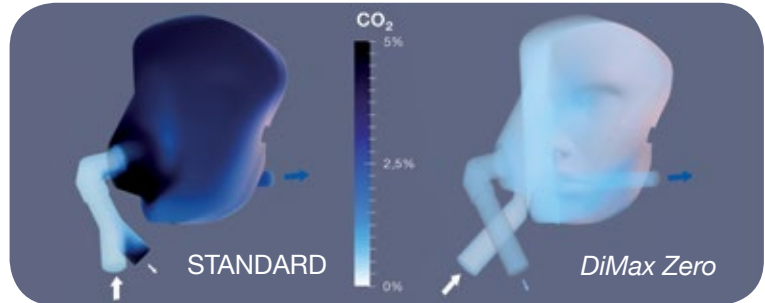

TECHNICAL FEATURES

- Zero dead space
- 5-point headgear attachment
- Larger nose room
- Dedicated ports for NG tube and bronchoscopy probe
- Single patient and reprocessable versions
- Interchangeable elbows:
 - 22F non-vented (blue)
 - 22M vented (clear)

ADVANTAGES

- Zero CO₂ rebreathing
- Improved patient/ventilator synchronism
- Improved airtight performance
- Improved comfort regardless of patient face contour or nose size
- Less headgear pressure
- Less risk of skin lesions
- NG tubes and bronchoscopy probes can be used without air leaks or disrupting the ventilation
- Single multi patient versions or reprocessable in autoclave @ 134° - Sterrad - ETO in autoclave @ 134° - gas plasma - ETO



The fifth attachment point allows a better seal on patients' forehead



Dedicated air tight NG tube port



“Non-Vented” elbow for single limb ventilator



Configuration for CPAP with “Easy Vent” Venturi Generator



Configuration for dual limb ventilator



Configuration for single limb ventilator or CPAP Generator

REFERENCES



CODES*		SIZE
SINGLE PATIENT**	REPROCESSABLE***	
700/6828	700/6824	SMALL
700/6829	700/6825	MEDIUM
700/6830	700/6826	LARGE
700/6831	700/6827	EXTRA LARGE

* Including “non-vented” blue elbow, “vented” clear elbow and 1 headgear
 ** single packed
 *** single packed

DIMAR

Via G. Galilei, 6 41036 Medolla -MO- I
 Tel: +39 0535 611336 Fax: +39 0535 611328
 Mail: dimar.srl@tiscali.it

