EasyID® PATIENT IDENTIFICATION



Positive Patient Identification Made Easy ~ Because We Care



TABLE OF CONTENTS

Foam Ba	nc	S
---------	----	---

Band On A Sheet (BOAS)

Recommendations

Polyester Wristbands

Alert Wristbands

Tyvek Wristbands

Thermal Wristbands

Paper Labels

Paper & Polyester Labels

Thermal Labels

Foam Identification Bands for Infants, Newborns & Preemies

Our complete line of soft foam bands offers a safe, comfortable and adjustable solution for growing babies with sensitive skin. They feature an adjustable Velcro® closure, ensuring baby's safety and comfort. The bands can be used in conjunction with individual printed thermal tags or with our Mother & Baby - Band On A Sheet (BOAS) solution. Our versatility allows facilities to attach different types of labels to help meet their needs for a reliable Positive Patient Identification. The bands with extensions make identifying a swaddled baby quick and easy by allowing staff to scan barcodes or read ID information outside of blankets. Our ID bands are free from latex and vinyl.



Product #: 20-321W Simple Infant Band



SHIELD: No SIZE: 1.25 x 5.5 inch COLOR: White

MATERIALS: Foam with Velcro®

closure

PACKAGING: 24 bands per box

Product #: 20-325W

Simple Pediatric Band



SHIELD: No SIZE: 1.25 x 7.75 inch COLOR: White

MATERIALS: Foam with Velcro® closure

closure

PACKAGING: 24 bands per box

Product #: 20-115W2

Thermal Tag Label



LAYOUT: Roll (1 inch inner core) SIZE: 1.25 x 2.25 inch per tag

COLOR: White

MATERIAL: Direct thermal tag AMOUNT: 500 labels per roll PACKAGING: 4 rolls per case

Product #: 20-185W3

Tag w/Band On A Sheet (BOAS)



LAYOUT: 8.5 x 11 inch sheet SIZE: 2.5 x 1.25 inch for tag

COLOR: White

MATERIAL: Polyester tag

AMOUNT: 1 tag per sheet

BANDS: 1 infant & 1 adult

LABELS: 20 - 1 x 2.5 inch

PACKAGING: 500 sheets per case

SIMPLE - Pediatric Foam Band #20-325W



SHIELDED - Infant Foam Bands #20-322 White - Blue - Pink



Product #: 20-322W Shielded Infant Band




SHIELD: Yes
SIZE: 1.25 x 5.5 inch
COLOR: White
MATERIALS: Foam with V

MATERIALS: Foam with Velcro® closure and shield

PACKAGING: 24 bands per box

Product #: 20-322B
Shielded Infant Band

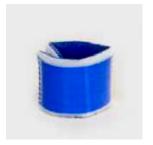

SHIELD: Yes

SIZE: 1.25 x 5.5 inch

COLOR: Blue

MATERIALS: Foam with Velcro® closure and shield

PACKAGING: 24 bands per box

Product #: 20-322P
Shielded Infant Band



SHIELD: Yes
SIZE: 1.25 x 5.5 inch
COLOR: Pink
MATERIALS: Foam with Velcro® closure and shield
PACKAGING: 24 bands per box

Product #: 20-328W Shielded Pediatric Band



SIZE: 1.25 x 7.75 inch
COLOR: White
MATERIALS: Foam with Velcro® closure and shield

PACKAGING: 24 bands per box



SHIELDED - Pediatric Foam Band #20-328W

Foam Identification Bands with Extensions

Our addition of specialized extensions allows for identifying a swaddled baby quickly and easily, without disrupting them from their sleep.

POLYESTER EXTENSION - Infant Foam Band #20-326W



Product #: 20-326W Poly Ext. Infant Band



SHIELD: No
SIZE: 1.25 x 5.5 inch
COLOR: White
MATERIALS: Foam, Polyes

MATERIALS: Foam, Polyester and Velcro® closure

PACKAGING: 24 bands per box

Product #: 20-323W Cloth Ext. Infant Band



SHIELD: No
SIZE: 1.25 x 5.5 inch
COLOR: White
MATERIALS: Foam, Cloth and

MATERIALS: Foam, Cloth and Velcro® closure

PACKAGING: 24 bands per box

Product #: 20-322WE
Cloth Ext. with Shielded Band



SHIELD: Yes

SIZE: 1.25 x 5.5 inch

COLOR: White

MATERIALS: Foam, Cloth, Velcro®

closure and clear plastic shield PACKAGING: 24 bands per box


Product #: 20-322WP
Poly Ext. with Shielded Band



SHIELD: Yes SIZE: 1.25 x 5.5 inch COLOR: White

MATERIALS: Foam, Poly, Velcro® closure and clear plastic shield PACKAGING: 24 bands per box



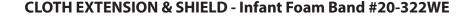



#20-115W2

Tag: #20-118W3 DOB: 01/01/20XX Age: 1D Sex: M DR. BARNES 231 MRN: 0007858 071423





Tag w/Band On A Sheet (BOAS)

Product #: 20-118W3

Poly Tag w/Paper Labels



LAYOUT: 8.5 x 11 sheet

SIZE: 1.875 x 1.3125 inch per tag

COLOR: White

MATERIAL: Polyester tag AMOUNT: 2 tags per sheet PACKAGING: 1,000 sheets per case

Product #: 20-115W2

Thermal Tag Label



LAYOUT: Roll (1 inch inner core)

SIZE: 1.25 x 2.25 inch per tag

COLOR: White

MATERIAL: Direct thermal tag AMOUNT: 500 labels per roll PACKAGING: 4 rolls per case

LAYOUT: 8.5 x 11 inch sheet SIZE: 2.5 x 1.25 inch per tag

COLOR: White

MATERIAL: Polyester tag AMOUNT: 1 tag per sheet BANDS: 1 infant & 1 adult LABELS: 20 - 1 x 2.5 inch

PACKAGING: 500 sheets per case

In addition to our adjustable Velcro® closure, we offer cloth and polyester extensions, along with a variety of waterproof tags, to eliminate staff from needing to unswaddle a baby during a barcode scanning to achieve a positive identification. The specially processed polyester

and thermal tags offer hospitals an easy solution to directly print onto tags and labels while protecting patient information from moisture and hand-sanitizers, preserving patient data for long hospital

stays and usability.

The band with a polyester extension can be used in conjunction with individual printed thermal tags or with our Mother & Baby - Band On A Sheet (BOAS) solution.

Our BOAS are free from latex, heavy plastics and vinyl.

Foam Band #20-322WE with BOAS: #20-185W3



EasyID® Band On A Sheet (BOAS)

Laser Sheets with Wristbands & Labels

EasyID®'s Laser Printed Band On A Sheet (BOAS) consists of polyester wristbands and labels that are made with a resilient paper material. The wristbands are by far the most comfortable, lightweight bands on the market, while the cost saving labels are easy to peel and adhere to patient charts and records. The BOAS integrates easily with existing laser printers, enabling hospitals to print wristbands and labels simultaneously from the same laser printer. All wristbands and labels are customizable and available in Adult, Pediatric, Infant, Mother and Infant, or Father and Infant sizes.





Product #: 20-181W3

Wristband on Bottom - Band On A Sheet (BOAS)

SHEET SIZE: 8.5 x 11 sheet BAND: 1 - 1 x 10.75 inch LABELS: 24 - 1 x 2.5 inch COLOR: White and Yellow

MATERIAL: Paper Labels and Polyester Wristband

* Replacement for LaserBand® PLS-102W and Standard Register® HLW20

PACKAGING: 500 sheets per case

MITH, HENRY J.

WAS DISTRICT AND 13 MIN 13 M

Product #: 20-189W3

Wristband on Top - Band On A Sheet (BOAS)

SHEET SIZE: 8.5 x 11 sheet BAND: 1 - 1 x 10.75 inch LABELS: 24 - 1 x 2.5 inch

COLOR: White

MATERIAL: Paper Labels and Polyester Wristband

PACKAGING: 500 sheets per case

ADULT BAND										
EasyI D ®										
		Labels								
	Chai									
			20-189W3							

EasyID® Band On A Sheet (BOAS)

By using the latest in ID band technology, we are able to provide a wristband on a sheet that doesn't require an extra plastic shield to protect printed text and/or images. Our BOAS are designed to replace LaserBand® and Standard Register® laser sheets.

EasyID®



Product #: 20-182W3*

Two Adult Wristbands - Band On A Sheet (BOAS)


SHEET SIZE: 8.5 x 11 sheet BANDS: 2 - 1 x 10.75 inch LABELS: 20 - 1 x 2.5 inch

COLOR: White

MATERIAL: Paper Labels and Polyester Wristbands

PACKAGING: 500 sheets per case

^{*} To order, please call for availability and requirements

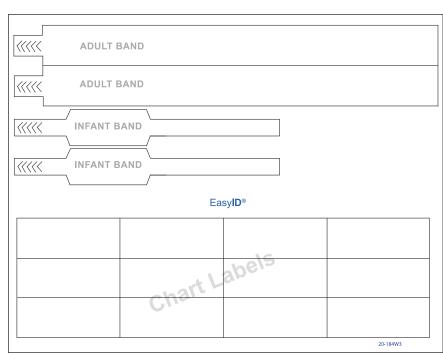



EasyID®

Product #: 20-184W3

Two Infant & Two Adult Wristbands- Band On A Sheet (BOAS)

SHEET SIZE: 8.5 x 11 sheet BANDS: 2 - 1 x 10.75 inch 2 - .625 x 6.6875 inch

LABELS: 12 - 1 x 2.625 inch

COLOR: White

MATERIAL: Paper Labels and Polyester Wristbands


PACKAGING: 500 sheets per case

⁰ **////** ADULT BAND \bigcirc ADULT BAND **/////**

EasyID® Band On A Sheet (BOAS)

Band On A Sheet (BOAS) for Infant & Adult

Offering the latest in laser wristband technology, our BOAS for infant and adult easily replaces other major brands and is standard when it comes to convenient solutions for patient security. The bands are waterproof, tamper evident, tear resistant and have an average life of up to 14 days. The material is excellent for barcode printing from laser printers, and is ideally suited for healthcare systems for patient care, laboratories, and admissions. Printing can be achieved on one single sheet that easily allows both wristbands and labels to be printed simultaneously from the same laser printing tray. The additional waterproof ID tag label can be used in conjunction with our foam bands (with and without extensions) for patients with extra sensitive skin, or can be used to meet required double banding for both baby and mother. All of these features allow for a one-stop process for a Positive Patient Identification.

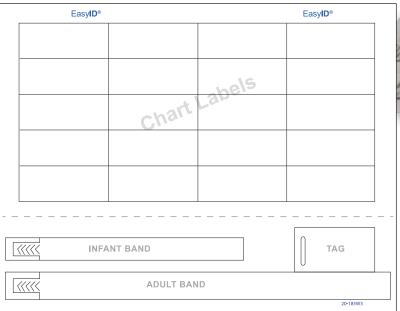



Product #: 20-185W3

Infant & Adult - Band On A Sheet

SHEET SIZE: 1-8.5 x 11 inch ADULT BAND: 1 - 0.75 x 10.5 inch INFANT BAND: 1 - 0.625 x 6.5 inch TAG LABEL: 1 - 2.5 x 1.25 inch

LABELS: 20 - 1 x 2.5 inch

COLOR: White

MATERIAL: Paper Labels and Polyester tag and bands

PACKAGING: 500 sheets per case

Laser printed ID band/label combination sheets are a popular option for many hospitals. As barcode scanning has increased, the need for a non-shielded ID band has evolved. Also, plastic shields reflect light, which hinders barcode scanning. Our BOAS feature a polyester ID band to eliminate the shield, which vastly improves machine readability of barcodes. Removing the shield also increases patient comfort and ease of use. Our sheets provide 20% more chart labels than our competitors, and the design of the BOAS also prevents jamming in printers.

KEY FEATURES & BENEFITS:

- Various designs
- Reduces skin irritation
- Soft, lightweight and durable
- Unshielded, no laminate required
- Replaces LaserBand®/Standard Register®
- Confident machine readability
- Extra chart labels
- 20% more room for patient ID
- Double adhesive closure
- Competitive pricing



EasyID® Recommendations

While our patient identification products are made to integrate seamlessly with most current hospital systems, there are many options to consider when selecting patient wristbands to meet your facility's needs. EasyID® Patient Identification Specialists will assist you to find the right solution. With so many wristband options and features to consider, it is important to evaluate your needs to achieve the most accurate patient ID solution.

Administration Department

To help streamline patient admissions, the 20-185W3 BOAS simultaneously prints patient chart labels and wristbands on one sheet from any standard laser printer. The sheet includes wristbands to fit all patients, from newborn to adult. The combination of labels and wristbands enables quick deployment of a barcoded solution without the need for separate printers and/or products. These easy-to-use wristbands do not require stick-on labels or inserts, which saves time and eliminates potential error.

Emergency (ER) | Short-Stay | Outpatient

Short-term hospital applications, such as the ER, require patient ID that is accurate, reliable, economical, and effective. Any of our BOAS provides a variety of wristbands and chart labels that can print from a standard laser printer to lessen the need for multiple printers, creating a more efficient workflow.

We also offer Tyvek $^{\circ}$ wristbands (20-306) and lightweight polyester wristbands (21-302) to be used with our $\frac{1}{2}$ sheet waterproof labels (20-105W). These two solutions can be easily customized for any clinical area.

NICU Area

Matching infant and mother ID wristbands are made possible with our 20-185W3 and 20-184W3 BOAS. These combination sheets of wristbands and labels provide a multiple banding solution that is often required by hospitals. The additional tag label on our 20-185W3 BOAS can be used in combination with our foam bands to enable caregivers to easily scan and verify the patient ID without the need to unswaddle, keeping patients more comfortable and less stressed.

If a patient is going directly to the NICU without the need to print multiple labels, a single waterproof tag label (20-115W2) can be printed and attached to a foam band (20-321W, 20-326W & 20-323W), allowing caregivers to verify the patient ID without the need to unswaddle or unnecessarily move the patient. Also, since one tag can be printed at a time, there will be less waste when a full sheet is not needed.

Renal/Kidney | Cancer | Heart Centers

Our softer, lightweight bands help prevent skin irritations while improving patient wellbeing by keeping them comfortable. Our specially-made polyester wristbands, with or without shields (21-203W & 21-302W), provide many possibilities for caregivers to correctly identify that patients are in the appropriate area. We also offer many label solutions to coordinate with the varying needs of each department.

Psychiatric Center

While a softer band material is ideal for most patients, this option often requires a much stronger and thicker band. We offer a shielded band (20-205) that measures 1.25 x 10 inches and is made of a much thicker polyester material. It also features a shield under which a label can be placed so it cannot be removed. This band is currently being used in several different psychiatric facilities throughout the country. The above band could be used with a single sheet of labels similar to the current layout of 20-101W2. This extremely durable band is almost impossible to remove, and the plastic shield allows the label to be permanently affixed to the band to prevent swapping of identification.

Some of our most popular label sheets to be used with shielded bands include the following:

20-090W2 20-091W2 20-101W2 20-103W 20-131W

EasyID® Polyester Wristbands

Polyester wristbands without a shield

Our unshielded polyester wristbands offer an excellent alternative to shielded wristbands. The bands are waterproof, tamper evident and tear resistant. By applying one of our polyester labels on top of these durable polyester bands, the labels are scannable by barcode scanners, allowing for an effective and efficient way to get a computer confirmation of the patient's identity and prevent hospital medical errors. Hospitals can use any label.

Product #: 20-310 Product #: 20-301

Pediatric Polyester Band without shield

Narrow Polyester Band without shield

Polyester Band without shield


Product #: 20-302 Product #: 21-302 Product #: 20-308

Polyester Band without shield

Long Polyester Band without shield




- .625 x 7 inch
- No shield
- 500 bands per box

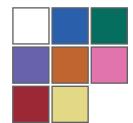

- .625 x 10 inch
- · No shield
- 500 bands per box



- 1 x 10 inch
- · No shield
- 500 bands per box

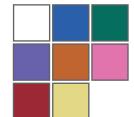


- 1 x 11 inch
 - · No shield
 - 5 bands on a sheet
 - 400 bands per box






- 1 x 13 inch
- No shield
- · 250 bands per box

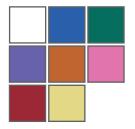


Polyester labels - ideal for wristbands without a shield



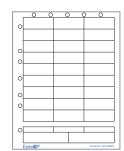


Product #: 20-109W3

27 Paper Chart Labels 5 Polyester Labels for unshielded wristbands

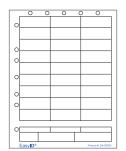

Product #: 20-120W3

27 Paper Chart Labels **6 Polyester Labels for** unshielded wristbands

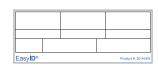

Product #: 20-105W

9 Polyester Labels for multiple use and unshielded wristbands

EasyID® Polyester Wristbands

Polyester wristbands with shields

Shielded polyester wristbands completely seal a paper ID label behind an impenetrable shield for the best protection against all kinds of liquids. This provides hospitals with a waterproof, high-performing, latex-free wristband for quality patient care and safety. These time-tested and proven performance bands, with a securely sealed adhesive closure for added safety, come in vibrant, solid colors for excellent visibility for additional identification.

Product #: 20-201

Product #: 20-303

Product #: 20-203

Product #: 21-203

Product #: 20-206

Pediatric Polyester Band with shield

Narrow Polyester Band with shield

Polyester Band with shield

Polyester Band with shield

Long Polyester Band with shield














- .625 x 7 inch
- · With shield
- 500 bands per box
- .625 x 10 inch
- · With shield
- 500 bands per box
- 1 x 10 inch
- · With shield
- 500 bands per box
- 1 x 11 inch
- · With shield
- 5 bands on a sheet
- · 400 bands per box

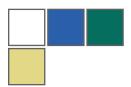





- · With shield
- 250 bands per box

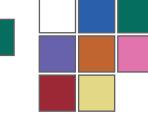

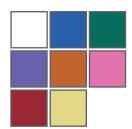


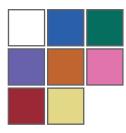

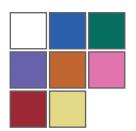

Product #: 20-205 Wide Polyester Band with shield



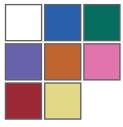



- 1.25 x 10 inch
- · With shield
- · 250 bands per box



EasyID® Alert Wristbands

Patient Alert Wristbands

In an effort to improve patient safety, the American Hospital Association (AHA) recommends all member hospitals apply a standardized color process of the following alert wristbands: red for allergy, yellow for fall risk and purple for do not resuscitate.

To further increase patient safety and comfort, EasyID® proudly introduces the 21-311 series. The new .625 x 11 inch band features a lightweight polyester material for maximum durability, including an extra strong adhesive for a secure closure.

Product #: 21-311RA3

Polyester Alert Band without shield



- .625 x 11 inch
- · No shield
- · Material: Polyester
- 400 bands per box

Product #: 21-311YF3

Polyester Alert Band without shield



- .625 x 11 inch
- No shield
- · Material: Polyester
- 400 bands per box

Product #: 21-311LD3

Polyester Alert Band without shield



- .625 x 11 inch
- No shield
- Material: Polyester
- 400 bands per box

Product #: 21-311GLA

Polyester Alert Band without shield



- .625 x 11 inch
- No shield
- · Material: Polyester
- 400 bands per box

Product #: 21-311PRL

Polyester Alert Band without shield



- .625 x 11 inch
- No shield
- · Material: Polyester
- 400 bands per box

Product #: 21-311LM

Polyester Alert Band without shield



- .625 x 11 inch
- No shield
- Material: Polyester
- 400 bands per box

Product #: 21-315SA

Polyester Alert Band without shield



- .625 x 17 inch
- No shield
- · Material: Polyester
- 400 bands per box

Product #: 20-311YF

Tyvek® Alert Band without shield



- .625 x 10 inch
- No shield
- Material: Tyvek®
- 800 bands per box



Call today for Custom Alert Bands

Tyvek® wristbands without a shield

Tyvek® ID bands are a great option for "short-stay" patient visits. They offer a lightweight band to ensure maximum patient comfort during stressful visits and are an excellent option for emergency and ambulatory departments. In addition, we provide a complete and effective positive identification solution with the use of our waterproof polyester labels.

Each wristband has a strong, waterproof adhesive closure and is made with a genuine DuPont Tyvek® material. They are a great economical choice for short-term use.

Product #: 20-304

Infant Tyvek® Band without shield



- .625 x 7 inch
- · No shield
- 200 bands per box



Product #: 20-312R

Narrow Tyvek® Band without shield



- .625 x 10 inch
- No shield
- 800 bands per box



Product #: 20-311YF

Tyvek® Alert Band without shield



- .625 x 10 inch
- · No shield
- 800 bands per box




Product #: 20-306


Tyvek® Band without shield



- 1 x 10 inch
- No shield
- 500 bands per box

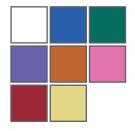


Polyester labels - ideal for wristbands without a shield

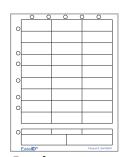

Product #: 20-109W3

27 Paper Chart Labels 5 Polyester Labels for unshielded wristbands

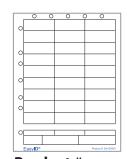

Product #: 20-120W3

27 Paper Chart Labels 6 Polyester Labels for unshielded wristbands



Product #: 20-105W

9 Polyester Labels for multiple use and unshielded wristbands



EasyID® Thermal Wristbands

Direct Thermal Wristbands - G series printer

Direct Thermal Latex-Free Wristbands are an effective, error-free patient identification solution that ensures patient safety with a point-of-care process. The wristbands have a strong adhesive backing, are waterproof and made of a tear-resistant synthetic material that allows hospitals to print anywhere. These thermal wristbands help hospitals prevent medical errors by printing accurate and reliable patient information.

Product #: Product #: Product #: Product #: Product #: 20-664W 20-663W 20-662W 20-661 20-651 **Infant Thermal Band Child Thermal Band Thermal Band Wide Thermal Band Narrow Thermal Band** without shield without shield without shield without shield without shield • 1 x 6 inch • 1 x 7 inch • .75 x 11 inch • 1 x 11 inch • 1.25 x 11 inch · No shield • 350 bands per roll • 300 bands per roll • 200 bands per roll · 200 bands per roll · 200 bands per roll 6 rolls per case 6 rolls per case 6 rolls per case 6 rolls per case Sold per roll

Printers for Direct Thermal Wristbands and Labels

Print your labels and/or wristbands on the go with our Thermal Printers:

We offer a complete line of stationary, desktop and mobile printers, along with accessories, customer support and supplies. Our printers are perfect for all of your light to medium duty thermal barcode printing applications and offer a great solution for printing on all types of specialty materials. They are engineered for ease-of-use while offering a broad range of connectivity options, making them ideal for countless low-to mid-volume label, tag or wristband printing applications.

Our series of lightweight, rugged mobile printers provides user-friendly options and flexible connectivity for any complex mobile printing application. They are built to handle increased levels of wireless security, making them ideal for secure wireless transfer of sensitive patient data which provides staff the ability to print patient identification onto labels and wristbands at point-of-care. Our EasyID® Mobile solutions leverage the power of these printers to provide an easy point-of-care patient ID confirmation and label printing solution.

Contact us today to determine the best thermal printer and ID system to meet your needs.

EasyID® Thermal Wristbands

Direct Thermal Wristbands - HC100 Cartridge

Thermal Wristbands can be printed with a linear and/or 2-D barcode that features high-resolution and high contrast print that can be scanned immediately. This high-resolution print, on a durable latex-free material, produces a precise printed barcode that delivers an immediate first-time scan while remaining smudge-proof and accurate throughout the patient's stay. Thermal wristbands are waterproof, soft, lightweight and smooth to the patient's skin. The unique design provides a level of comfort and is latex-free, hypoallergenic, and does not contain skin-irritating antimicrobial agents or other harsh chemical additives that can cause allergic reactions to patients.

Product #:

Thermal Band

without shield

20-680

Product #: 20-683WC

Soft Infant Band without shield



- 1 x 7.6875 inch
- · No shield
- 275 bands per cartridge
- 6 cartridges per case



Product #: 20-682WC

Comfort Infant Band without shield



- 1 x 6 inch
- · No shield
- 600 bands per cartridge
- 6 cartridges per case



- 1 x 11 inch
 - · No shield
 - 200 bands per cartridge
 - 6 cartridges per case



Product #: 20-681WC

Comfort Pediatric Band without shield



- 1 x 7 inch
- · No shield
- 500 bands per cartridge
- 6 cartridges per case



Product #:

20-680WC

without shield

Comfort Adult Band

- 1 x 11 inch · No shield
- 300 bands per cartridge
- · 6 cartridges per case







Direct Thermal Wristbands are printed from a small portable thermal printer, which allows you to print patient information onto a band instantly. There are no extra labels to print and apply to the wristbands, ensuring a one-step process. Within seconds, you can print and apply the wristband on a patient. Furthermore, the small space saving printer is easy to install, which eliminates the need to destroy ribbons that contain important patient information.


Call today to receive a quote on these small portable printers.

Paper Labels with Multiple Configurations

Hospitals are continuously under pressure to save time and money without sacrificing accuracy or quality. We help hospitals evaluate their daily operations as it relates to applying label systems and products to ensure patient safety. Single sheet printing provides hospitals with an easy one-step label and wristband solution that can save time and money. Paper Label Sheets and Bands On A Sheet (BOAS) are the most widely used labels because of their low cost and ability to print using a variety of printers, systems and materials.

Product #: 20-090W2

Patient paper labels on standard size sheet

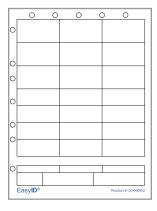

SHEET SIZE: 8.5 x 11 inch sheet

HOLES: 12 - 5 on top, 7 on side

LABELS:

9 - 1.375 x 2.5 inch paper labels

12 - 1 x 2.5 inch paper labels

3 - 0.4375 x 2.5 inch paper labels

2 - 0.75 x 3 inch paper labels

1 - 0.75 x 1.5 inch paper label

MATERIAL: Paper


PRINTER: Standard Laser

PACKAGING: 1,000 sheets per box

COLOR: White

Product #: 20-150W2

Patient paper labels on standard size sheet

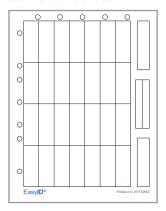

SHEET SIZE: 8.5 x 11 inch sheet

HOLES: 12 - 5 on top, 7 on side

LABELS:

28 - 0.94 x 2.55 inch paper labels

2 - 0.4375 x 3 inch paper labels

2 - 0.75 x 3 inch paper labels

MATERIAL: Paper


PRINTER: Standard Laser

PACKAGING: 1,000 sheets per box

COLOR: White

Product #: 20-179W2

Patient paper labels on standard size sheet

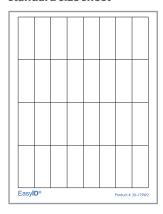

SHEET SIZE: 8.5 x 11 inch sheet

HOLES: None

LABELS:

32 - 1 x 2.625 inch paper labels

MATERIAL: Paper


PRINTER: Standard Laser

PACKAGING: 1,000 sheets per box

COLOR: White

Product #: 20-093W2

Patient paper labels on standard size sheet

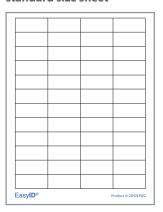

SHEET SIZE: 8.5 x 11 inch sheet

HOLES: None

LABELS:

48 - 2 x 0.875 inch paper labels

MATERIAL: Paper


PRINTER: Standard Laser

PACKAGING: 1,000 sheets per box

COLOR: White

Product #: 20-131W

Patient paper labels on 1/2 size sheet - A5 size

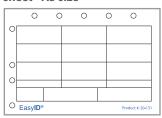

SHEET SIZE: 5.8 x 8.25 inch sheet 9 - 1 x 2.5 inch paper labels

HOLES: 9 - 5 on top, 4 on side

PRINTER: Standard Laser

MATERIAL: Paper

3 - 0.4375 x 2.5 inch paper labels

2 - 0.75 x 3 inch paper labels

1 - 0.75 x 1.5 inch paper label

PACKAGING: 1,000 sheets per box


COLOR: White

Paper Labels with Multiple Configurations

Paper labels are available in several different design configurations and in three sizes – full sheet, half sheet, and A5 sheet. Paper Labels can also be used as asset tags with printed barcodes on them. These can then be affixed to supplies and instruments in the hospital or medical clinic to provide positive identification and tracking, which assists in the process of checking hospital inventory supplies. EasyID® offers a comprehensive product line designed for use with any positive patient ID/forms automation system.

Product #: 20-092W2

Patient paper labels on standard size sheet

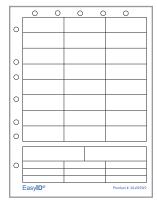

SHEET SIZE: 8.5 x 11 inch sheet

HOLES: 12 - 5 on top, 7 on side

LABELS:

24 - 0.95 x 2.5 inch paper labels

9 - 0.4375 x 2.5 inch paper labels

2 - 0.875 x 3.75 inch paper labels

MATERIAL: Paper

PRINTER: Standard Laser


PACKAGING: 1,000 sheets per box

COLOR: White

*To order, please call for availability and requirements

Product #: 20-151W2

Patient paper labels on standard size sheet

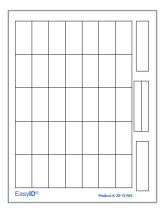

SHEET SIZE: 8.5 x 11 inch sheet

HOLES: None

LABELS:

35 - 2 x 1 inch paper labels

2 - 0.4375 x 3 inch paper labels

2 - 0.75 x 3 inch paper labels

MATERIAL: Paper

PRINTER: Standard Laser

PACKAGING: 1,000 sheets per box


COLOR: White

*To order, please call for availability

and requirements

Product #: 20-126W2

Patient paper labels on standard size sheet

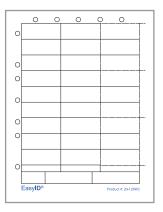

SHEET SIZE: 8.5 x 11 inch sheet

HOLES: 12 - 5 on top, 7 on side

LABELS:

27 - 1 x 2.5 inch paper labels

3 - 0.4375 x 2.5 inch paper labels

2 - 0.75 x 3 inch paper labels

1 - 0.75 x 1.5 inch paper label

Perforations allow for:

3 sets of 3 -1 x 2 ½ inch paper labels

MATERIAL: Paper

PRINTER: Standard Laser

PACKAGING: 1,000 sheets per box

COLOR: White

*To order, please call for availability and requirements

Product #: 20-103W

Patient paper labels on 1/2 size



SHEET SIZE: 8.5 x 5.5 inch sheet

MATERIAL: Paper

PRINTER: Standard Laser

9 - 1 x 2.5 inch paper labels

2 - 0.75 x 3 inch paper labels

1 - 0.75 x 1.5 inch paper label

YOUR CUSTOM LABEL SHEET **HERE**

If you don't see the design and layout you need, call and let us design one for you. It's easy!

sheet


PACKAGING: 1,000 sheets per box **HOLES: None** 3 - 0.4375 x 2.5 inch paper labels COLOR: White

Paper Labels with Multiple Configurations

The selective use of color on patient labeling can be a valuable resource for providing patient safety. When used appropriately, color can enhance processes to add clarity to departments, medication dispensing and administration operations. Our labels are printed with fade resistant ink to prevent premature fading. We offer several label sheets with different color layouts to help reduce the likelihood of human error and quickly alert staff and physicians to conditions and warnings of patients.

Product #: 20-091W2

Patient paper labels on standard size sheet

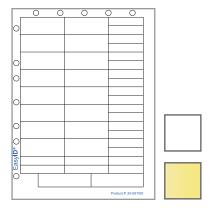

SHEET SIZE: 8.5 x 11 inch sheet

HOLES: 12 - 5 on top, 7 on side

LABELS:

16 - 1 x 2.5 inch paper labels

18 - 0.5 x 2 inch paper labels

4 - 0.5 x 2.5 inch paper labels

2 - 0.75 x 3 inch paper labels

MATERIAL: Paper


PRINTER: Standard Laser

PACKAGING: 1,000 sheets per box

COLOR: White and Yellow

Product #: 20-101W2

Patient paper labels on standard size sheet

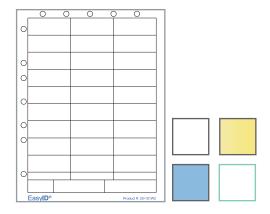

SHEET SIZE: 8.5 x 11 inch sheet

HOLES: 12 - 5 on top, 7 on side

LABELS:

27 - 1 x 2.5 inch paper labels

3 - 0.4375 x 2.5 inch paper labels

2 - 0.75 x 3 inch paper labels

1 - 0.75 x 1.5 inch paper label

MATERIAL: Paper

PRINTER: Standard Laser

PACKAGING: 1,000 sheets per box


COLOR: White, Blue, Yellow and

Green Border

According to a report by the Institute of Medicine, medical error is the eighth leading cause of death in patients in the United States. Improper labeling has caused many of these deaths due to overdoses, unnecessary surgeries and incorrect drugs being administered. EasyID® is reaching out to hospitals and medical facilities to help them improve Positive Patient ID (PPID) labeling by developing and customizing Patient Identification Systems, Labels, Forms, and Wristbands that will work together to improve workflow and patient safety. Building a reliable and easy Patient Identification Label System requires the right tools that will work with large and small hospitals and their medical ADT systems. EasyID® is making it possible for hospitals to easily implement a dependable identification system no matter where they are in the evolution of their patient information system.

Product #: 20-102W2

Patient paper labels on standard size sheet

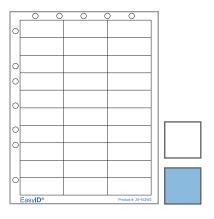

SHEET SIZE: 8.5 x 11 inch sheet

HOLES: 12 - 5 on top, 7 on side

LABELS:

30 - 1 x 2.5 inch paper labels

MATERIAL: Paper

PRINTER: Standard Laser

PACKAGING: 1,000 sheets per box

COLOR: White and Blue




Paper and Polyester Labels with Multiple Configurations

The combination of polyester and paper labels on the same sheet offers many cost effective solutions. While the paper labels are used for charts and documents, the polyester labels consist of a polymer material that hospitals can use to make accurate, durable wristbands. Our polyester labels are resistant to scuffing, tearing, abrasion, high heat, low temps, sunlight, chemical solvents, water and controlled shrinkage and can last a long time. Our paper and polyester label sheets are usually white, although other custom colors are available.

Product #: 20-109W3

Patient paper and polyester labels on standard size sheet

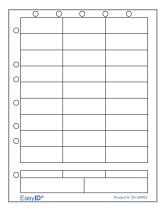

SHEET SIZE: 8.5 x 11 inch sheet

HOLES: 12 - 5 on top, 7 on side

LABELS:

27 - 0.95 x 2.5 inch paper labels

3 - 0.4375 x 2.5 inch poly labels

2 - 0.875 x 3.75 inch poly labels

MATERIAL: Paper and polyester

PRINTER: Standard Laser

PACKAGING: 1,000 sheets per


box

COLOR: White and Yellow



Product #: 20-120W3

Patient paper and polyester labels on standard size sheet

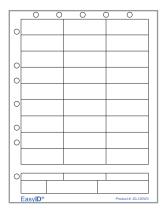

SHEET SIZE: 8.5 x 11 inch sheet

HOLES: 12 - 5 on top, 7 on side

LABELS:

27 - 0.95 x 2.5 inch paper labels

3 - 0.4375 x 2.5 inch poly labels

2 - 0.75 x 3 inch poly labels

1 - 0.75 x 1.5 inch poly label

MATERIAL: Paper and polyester

PRINTER: Standard Laser


PACKAGING: 1,000 sheets per

box

COLOR: White

Product #: 20-121W3

Patient paper and polyester labels on standard size sheet

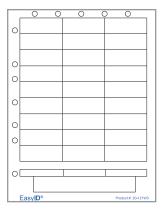

SHEET SIZE: 8.5 x 11 inch sheet

HOLES: 12 - 5 on top, 7 on side

LABELS:

27 - 0.95 x 2.5 inch paper labels

3 - 0.4375 x 2.5 inch poly labels

1 - 0.875 x 6 inch poly label

MATERIAL: Paper and polyester

PRINTER: Standard Laser


PACKAGING: 1,000 sheets per

box

COLOR: White

Product #: 20-118W3

Patient paper and polyester labels on standard size sheet

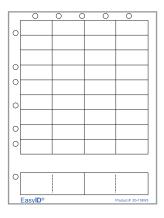

SHEET SIZE: 8.5 x 11 inch sheet

HOLES: 12 - 5 on top, 7 on side

LABELS:

36 - 0.875 x 1.875 inch paper labels

2 - 1.25 x 3.75 inch poly labels

MATERIAL: Paper


PRINTER: Standard Laser

PACKAGING: 1,000 sheets per box

COLOR: White

Product #: 20-105W

Patient polyester labels on 1/3 size sheet

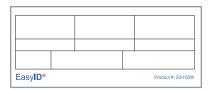

SHEET SIZE: 8.5 x 4 inch sheet

HOLES: None

MATERIAL: Polyester

PRINTER: Standard office

3 - 1 x 2.5 inch poly labels

3 - 0.5 x 2.5 inch poly labels

2 - 0.75 x 3 inch poly labels

1 - 0.75 x 1.5 inch poly label

PACKAGING: 1,000 sheets per box

COLOR: White

EasyID[®] Labels

Thermal Labels on a Roll



Product #	W (inch)	L (inch)	Inner Core (inch)	Outer Dia. (inch)	Printer Type	Matl. Type	Description	Adhesive Type	Perforated	Color	Labels per roll	Rolls per case
20-500W	0.5	2	1	5	Direct Thermal Desktop	Paper	Standard, coated direct thermal paper label with an all-temp permanent acrylic adhesive. Provides the optimal balance between performance and price. Coated with an all-temp permanent adhesive.	All-Temp Permanent Adhesive	Yes	White	2,500	12
20-463W	0.625	2.5	1	4	Direct Thermal Printer	Paper	Standard, coated direct thermal paper label with an all-temp permanent acrylic adhesive. Provides the optimal balance between performance and price. Coated with an all-temp permanent adhesive.	All-Temp Permanent Adhesive	Yes	White	2,000	1
20-461W	2	0.75	1	4	Direct Thermal Printer	Paper	Standard, coated direct thermal paper label with an all-temp permanent acrylic adhesive. Provides the optimal balance between performance and price. Coated with an all-temp permanent adhesive.	All-Temp Permanent Adhesive	Yes	White	1,750	6
20-117W	2	1	3	8	Direct Thermal Tabletop	Paper	Standard, coated direct thermal paper label with an all-temp permanent acrylic adhesive. Provides the optimal balance between performance and price. Recommended for indoor use.	All-Temp Adhesive	Yes	White	5,500	8
20-701W	2	1	1	5	Thermal Transfer Desktop	Poly	Premium white gloss thermal transfer polyester label with a high-performance permanent acrylic adhesive. Ideal for applications requiring up to 3 years outdoors, resistance to moderate chemicals and temperature exposure up to 300°F.	Permanent Adhesive	Yes	White	2,530	8
20-440W	2	1	1	5	Direct Thermal Desktop	Paper	Standard, coated direct thermal paper label with an all-temp permanent acrylic adhesive. Provides the optimal balance between performance and price. Recommended for indoor use.	All-Temp Adhesive	Yes	White	2,340	6
20-114W	2	1.25	0.75	2	Direct Thermal Mobile	Poly	Direct thermal matte polypropylene label with an all-temp permanent acrylic adhesive. Offers good resistance to scratching, smearing and moisture. Provides some flexibility for labeling curved surfaces.	All-Temp Adhesive	No	White	280	36
20-470W	2.25	0.4375	1	3	Direct Thermal Printer	Paper	Standard, coated direct thermal paper label with an all-temp permanent acrylic adhesive. Provides the optimal balance between performance and price. Coated with an all-temp permanent adhesive.	All-Temp Permanent Adhesive	Yes	White	500	1
20-445W	2.25	0.75	1	5	Direct Thermal Desktop	Paper	Direct thermal paper label with an all-temp permanent acrylic adhesive. Specially coated to provide optimal performance. Recommended for indoor use.	Permanent Adhesive	Yes	White	3,315	12
20-447W	2.25	1.25	1	5	Direct Thermal Desktop	Poly	Direct thermal matte polypropylene label with an all-temp permanent acrylic adhesive. Offers good resistance to scratching, smearing and moisture. Provides some flexibility for labeling curved surfaces.	All-Temp Adhesive	Yes	White	2,000	6
20-446W	2.25	1.25	1	5	Direct Thermal Desktop	Paper	Premium, coated direct thermal paper label with an all-temp permanent acrylic adhesive. Specially coated to provide optimal performance. Recommended for indoor use.	Permanent Adhesive	Yes	White	2,100	12
20-441W	3	1	1	5	Direct Thermal Desktop	Paper	Premium, coated direct thermal paper label with an all-temp permanent acrylic adhesive. Specially coated to provide optimal performance. Recommended for indoor use.	All-Temp Adhesive	Yes	White	2,340	6

EasyID® Labels

Product #	W (inch)	L (inch)	Inner Core (inch)	Outer Dia. (inch)	Printer Type	Matl. Type	Description	Adhesive Type	Perforated	Color	Labels per roll	Rolls per case
20-460W	3.5	1	3.25	7.4905	Direct Thermal Printer	Paper	Standard, coated direct thermal paper label with an all-temp permanent acrylic adhesive. Provides the optimal balance between performance and price. Coated with an all-temp permanent adhesive.	All-Temp Permanent Adhesive	Yes	White	4,540	4
20-450W	3.75	0.875	1	5	Direct Thermal Printer	Paper	Direct thermal synthetic paper label with an all-temp permanent acrylic adhesive. Provides the optimal balance between performance and price. Coated with an all-temp permanent adhesive.	All-Temp Permanent Adhesive	Yes	White	1,000	1
20-443W	4	1.5	1	5	Direct Thermal Desktop	Paper	Premium, coated direct thermal paper label with an all-temp permanent acrylic adhesive. Specially coated to provide optimal performance. Recommended for indoor use.	All-Temp Adhesive	Yes	White	1,620	6
20-700W	4	2	1	5	Thermal Transfer Desktop	Paper	Premium thermal transfer paper label with a permanent acrylic adhesive. Specially coated to provide optimal performance. Recommended for indoor use.	Permanent Adhesive	Yes	White	1,370	6
20-442W	4	2.5	1	5	Direct Thermal Desktop	Paper	Premium, coated direct thermal paper label with an all-temp permanent acrylic adhesive. Specially coated to provide optimal performance. Recommended for indoor use.	All-Temp Adhesive	Yes	White	1,000	6
20-702W	4	4	3	8	Thermal Transfer Tabletop	Poly	Thermal transfer polypropylene label with an all-temp acrylic adhesive for blood bag and IV bag labeling. Compliant with FDA 175.105	All-Temp Adhesive	Yes	White	1,110	4

EasyID® provides Thermal Labels as both Direct Thermal Labels (image created by label media) and Thermal Transfer Labels (image created by transfer of print from a thermal ribbon) to work with related thermal printers. They are available in many different sizes and are compatible with table top, desk top, mobile and portable printers. Please contact us for more information on additional sizes and custom products.

DIRECT THERMAL LABELS AND PRINTERS:

Direct thermal printers require the use of heat activated thermal paper and/or polyester labels, but require no transfers of ribbon. The print color is limited to black. Direct thermal printing is popular when the life of a label is less than 1 year. The primary benefit is an overall lower cost and easier media changeover and increased security since a ribbon is not required.

- · No need to change ribbons.
- · No need to discard ribbons for patient security.
- · Labels print only in black, unless they are preprinted.
- Labels have a shelf life of 1 year or less.
- · Labels are subject to heat or light.
- Labels are subject to abuse or friction.

THERMAL TRANSFER LABELS AND PRINTERS:

Thermal transfer printers require a carbon ribbon, which gets transferred onto a label by heat, hence the name "thermal transfer". The printing is very crisp and durable so this method is excellent for higher density barcodes and labels that require dependable scanning for patient information. Additional benefits of Thermal Transfer printing include the ability to print on paper or polyester surfaces.

- Labels can be printed in colors other than black.
- Labels have a shelf-life of over 1 year.
- · Labels are subject to heat or light.
- Labels are subject to abuse or friction.
- · Labels are printed with higher density barcodes.

KEEP EVERYONE HAPPY

WITH RELIABLE AND EASY-TO-USE PRODUCTS

EASY

The innovative design and customization of our products help with the overall efficiency of printing, applying and scanning each of our products. The EasyID® products provide easy-to-scan barcodes – reducing frustration and time spent re-scanning barcodes.

SAFE

Our products are latex free and free from sharp edges that can irritate sensitive skin. Furthermore, our easy-toscan barcodes will increase patient safety by producing accurate patient information.

COMFORTABLE

EasyID® products are comfortable and lightweight – you'll forget you are wearing a wristband. As a result, the quality of care is improved and people are happier.



LET US CREATE THE PERFECT PATIENT IDENTIFICATION SOLUTION FOR YOU TODAY!



7101 Executive Center Drive, Suite 123 Brentwood, TN 37027

866.371.5132

www.EasyID.org | www.gbscorp.com