

Olympic Papoose Board Accessories

Head immobilizer:

Complete immobilization of head. Attaches quickly to board. Adjustable Velcro straps. Foam pads for comfort. When ordering, specify Papoose Board size.



Flap extenders:

For immobilizing larger children. Two wide straps, with Velcro adjustment, attach to chest and abdominal flaps of Regular Size Papoose Board. Permits greater adjustment range.



Arm immobilizer:

For Small and Regular Size Papoose Boards. Slides under board to firmly secure arms. Velcro arm straps adjust laterally.



All Olympic Papoose Boards are latex-free.
Olympic Papoose Boards are fully guaranteed for one year.

Ordering information

Olympic Papoose Boards

Small Size For infants up to 24 months
Includes head & arm straps
CAT. NO. 50510

Regular Size For children ages 2 to 6 years
Includes head & arm straps
CAT. NO. 50500

Large Size For children 6 to 12 years
Includes head & arm straps
CAT. NO. 53600

Extra-Large Size For teenagers & adults
Includes head & arm straps
CAT. NO. 54610

Head Immobilizers

For Small Size Papoose Board
CAT. NO. 50518

For Regular Size Papoose Board
CAT. NO. 50508

For Large Size Papoose Board
CAT. NO. 53608

Other Accessories

Arm Immobilizers
For Regular & Small Size Papoose Boards
CAT. NO. 50509

Flap Extenders
For Regular Size Papoose Board
CAT. NO. 50524

Safe and secure immobilization –
Ages 3 months to teenagers and adults



Used in thousands of clinics, hospitals and dental offices

For safe immobilization of agitated patients



Olympic Papoose Boards are indispensable aids in thousands of hospitals and in doctors' & dental offices

Quick & easy to apply

- Three sets of opposing canvas flaps with Velcro® fasteners fold over the patient's body
- The diagonal shape of each board causes the flaps to criss-cross, providing a snug fit

Easy to clean

- Simply wipe the board with liquid disinfectant or soap and water
- Flap sets and arm straps can be machine washed

Selective immobilization permits exposure while maintaining control

- Any of the pairs of flaps may be opened while the child is still held securely by remaining flaps and/or head & arm straps



Hangs on wall For compact storage.

Infinite adjustment Generous areas of Velcro on flaps.

Detachable flaps Easily cleaned.

Rigid Board Lightweight, durable. Wipes clean.




Head Strap Adjusts up or down. Foam padded.

Arm Holes Permit exposure of one or both arms.

Arm Straps With secure Velcro closure. Adjust up or down.

Diagonal Shape Makes flaps conform to body shape.

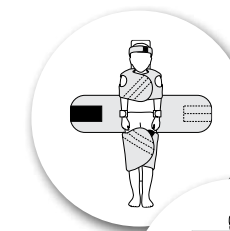


Abdomen Exposed

Chest and leg flaps closed. Arms and head held by straps.




Upper Body Exposed

Leg flaps closed. Arms immobilized by Velcro arm straps.

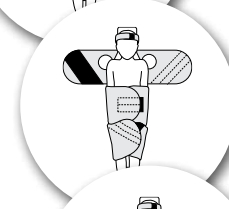


Chest Exposed

Abdominal and leg flaps closed to securely hold patient.

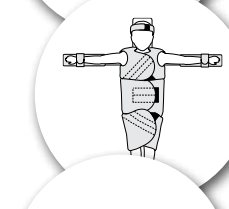


Arms Exposed

Through opening in chest flaps. Can be held on arm immobilizer (optional).

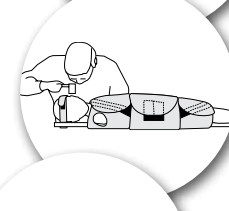


Head Strap

Adjustable Velcro. Can move up or down. (Optional head immobilizer available.)

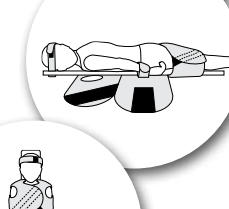


Prone Position

Leg flaps closed. Arms held by straps. Head secured with head strap.

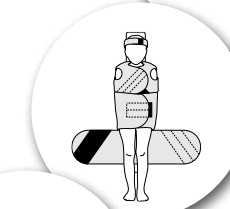

Legs Exposed

Chest and abdominal flaps closed for body immobilization.

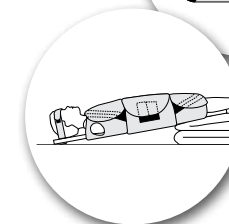

Trendelenburg Position or Reverse Trendelenburg.