



COMPETITORS






- Transparent needle free port (high rate);
- Transparent spike;
- 0,22-micron filter (alcohol-fobic);
- Smart skirt 20 and 13 mm (removable)
- Cap





SPIKES



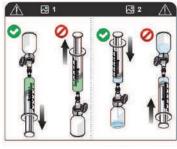









- Anti-disconnection system;
- Spinning 360°;
- Compatibility 100% needle free;
- High flow rete;
- Minimum dead space;
- Protected cap

COMPETITORS





















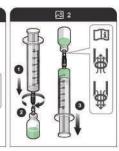

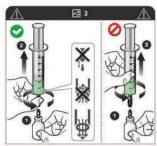

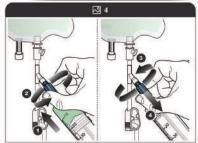

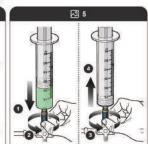



Standard sets

Our standard systems eliminate the risk of drugs spreading during transport and protects the operator from contamination.



- Break off sets
- A patent-protected system that guarantees the closure of the circuit until it is activated by the operator after connection to the administration set. This system eliminates the risks of drug spreading during transport and protects the operator.





- Y-needle free port (high rate);
- Amber PVC DHEP free tubing;
- Red Clamp (hazardous drug);
- Rotating luer lock with filter (filling line)

Buddy sets

Patent-protected anti-disconnection system that allows irreversible connection to the administration set, preventing subsequent disconnection from it.




- Amber PVC DHEP free tubing;
- Blu clamp (not hazardous drug);
- 2/4 one-way valves to prevent the reflux of chemotherapy;
- Adaptor for pump IV set or IV set.







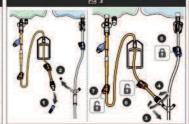





ADMINISTRATION


COMPOUNDING

- Transparent needle free port (high rate);
- 0,22 micron filter;
- Cap;
- Transparent needle free port (high rate);
- Protected needle;

BIOMEDICAL SOLUTIONS

- Easy insertions to secondary part of the bag;
- Allow to the European directive 2010/32/UE



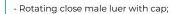


- Amber PVC DHEP free tubing;;
- Adaptor for pump IV set or IV set.










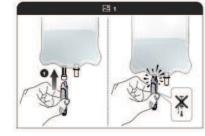

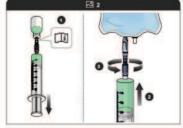








Click Lock™ System

CLICK LOCK™ SYSTEM is a new hybrid system which allows you to work through standard connections (needle-free luer) but also with an easy connection, guaranteed by a click.


When our **close male luer** comes in contact with the **neutral pressure valve** the connection is fast and avoids wrist rotation, allowing operators to work quickly and safely.




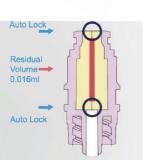





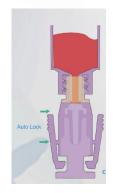

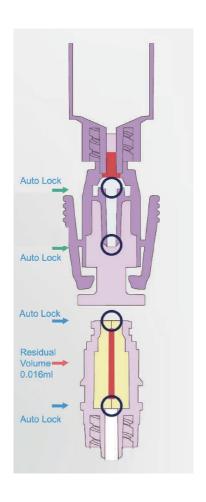

CLICK LOCK™ SYSTEM vs luer lock:

- Safer
- Easier to use
- More ergonomic
- 100% swabable
- Audible indicators
- No drip
- Hybrid

CLICK LOCK™ SYSTEM vs CSTD:

- Cheaper
- Easier to use
- Customizable